

Change Notice No.	Notification date	Implementation date
CN/19/004	June 05, 2019	July 05, 2019

Subject : New design of vent/drain valves of mdi Large Capsule Filters

Scope : This change notification will affect the products with following catalog numbers: (# can be any alphabet or numeral):

- LB#####
- LD#####
- LG#####
- LK#####
- LN#####
- LO#####
- LP#####
- LT#####
- LV#####
- LW#####

Background:

In our efforts to provide the user better products which are more convenient to use and deliver better performance, certain changes have been made in vent/drain valves of **mdi** Large Capsule Filters with catalog numbers as mentioned above.


Need for Change:

The existing vent/drain valves of **mdi** Large Capsule Filters with catalog numbers as mentioned above are validated to perform the venting operation during use. After venting operation, these vent/drain valves need to be closed/tightened in order to prevent leakage through vent/drain valves. It was observed that when existing vent/drain valves were closed using excessive force, it sometimes used to get deformed/damaged resulting into leakage through the same. To overcome this problem/issue, the vent/drain valves of **mdi** Large Capsule Filters have been made more robust through improvement in design. The new design of vent/drain valves has been verified for critical product performance parameters including functionality and ability to withstand excessive force during venting operation. No change has been done in any other critical dimensions or materials of construction or manufacturing process of **mdi** Large Capsule Filters.

Change Notice No.	Notification date	Implementation date
CN/19/004	June 05, 2019	July 05, 2019

How Does It Affect the User:

From a practical point of view the changes done will improve the ability of vent/drain valves to withstand excessive force during venting operation of **mdi** Large Capsule Filters. All performance specifications and product attributes will remain unchanged, except the improved design of vent/drain valves.

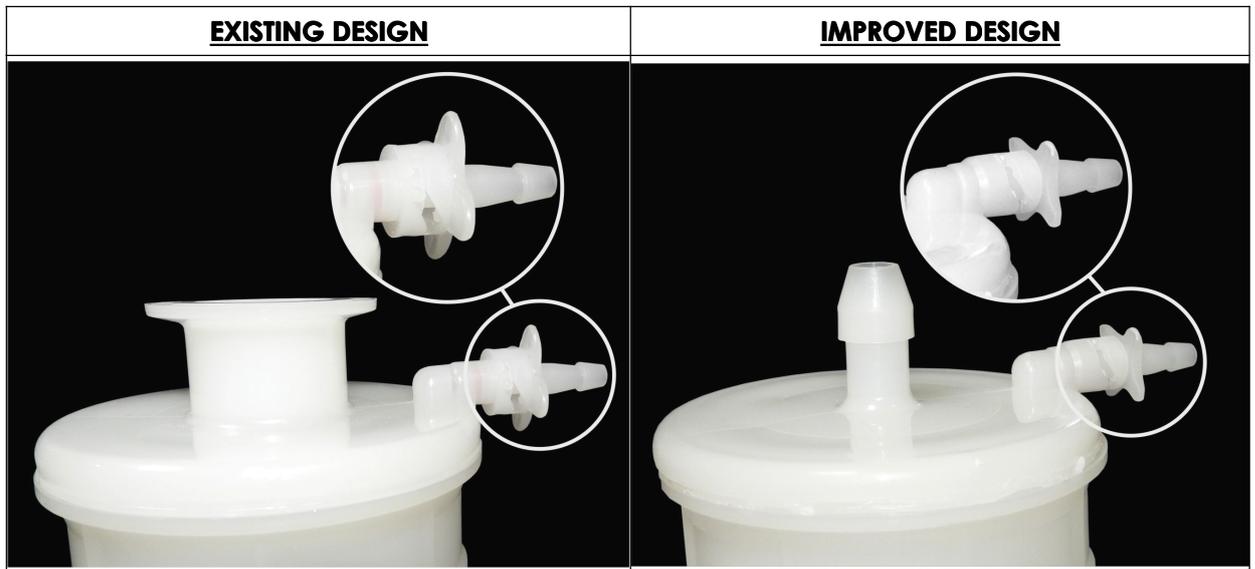

Implementation of Change:

mdi Large Capsule Filters having catalog numbers as mentioned above, manufactured after the implementation date (**July 05, 2019**) will have the said changes. However, the available stock of **mdi** Large Capsule Filters with catalog numbers as mentioned above having vent/drain valves with old design will be received by you till the stock lasts.

In case you have any queries, please feel free to contact our Technical Support Team at 'info@mdimembrane.com'.

Head - Quality Assurance
Advanced Microdevices Pvt. Ltd.