




# Easy Printer

Model LPM-02

The dispensing of reagents on the membranes has never been so simple and cost effective. The *Easy Printer* (Model : LPM-02) is the printer of choice for laboratory development of rapid assays. It is convenient not only for printing occasional lines with minimal reagent loss but due to in-built automation, can also be used for small pilot production.

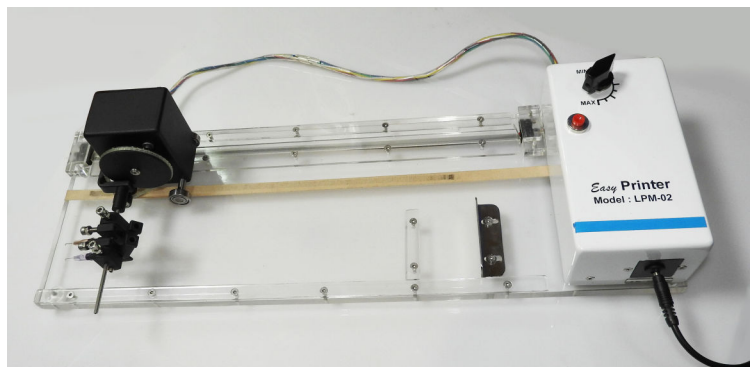

## SPECIFICATIONS

Maximum Print Length	30cm
No. of Printing Heads	2
Minimum Working Volume	50 µl
Voltage	220V, 50Hz/ 110V, 60Hz
Power	10W
Weight	3 Kg
Dimension	51 x 16.5 x 11.5 cms.
Dimension: Shipping Carton	51 x 21 x 20 cms.

## FEATURES

- Uniform dispensing of reagents results in good line quality.
- Two lines can be printed simultaneously
- Provision of variable speed for controlled dispensing of reagents.
- Auto-reverse facility makes the operation automatic.
- Pressure on the membrane is adjustable by a screw in the tip holder.
- Upto 600µl in holding tank allows printing of 20 laminates before refilling.
- Total recovery of unused reagents.

## ORDERING INFORMATION

*Easy Printer* Model LPM-02. Please indicate the operating voltage.

## ADVANCED SENSOR SYSTEMS PVT. LTD.

20, INDUSTRIAL AREA, AMBALA CANTT – 133 006, INDIA  
TEL: +91-171-2699471, Website: [www.mdmembrane.com](http://www.mdmembrane.com)  
Email: [advancedsensorysystemsambala@gmail.com](mailto:advancedsensorysystemsambala@gmail.com)